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1. Please fully charge before you use it. (Unplug the power before use, otherwise it won't open)
2. Massage head is faced down, and presses the button for few second to turn on the related functions.
Button H  : FDA-reg heat therapy, photobiomodulation
Button D  : Deep tissue stimulation, photobiomodulation
Caution: Do not directly view the light.
3. Press on the painful area gently and horizontally. Gently circle around it and do not stay at the same point for a long time.

It is recommended to apply the medical ointment after the use of MT Alpha, so that the medicine can be absorbed more effectively with the help of warm effect. Suggested to use MTcure's medicament.
4. Treatment level
Button H  : After switch on, the middle mode is default setting, press it to become high mode.
Button D  : After switch on, the middle mode is default setting, press it to become high mode, press it again to become light mode.
Caution: Do not directly view the light.
5. Please wipe the massage head with a cloth dampened with water or alcohol after use. Suggested to use 15-20 minutes and 1-2 times every day. Automatically turned off after 10 minutes.

Remarks:


- When the power is about to run out, the charging indicator will flash to remind the user to charge.
- It is recommended to switch off the machine after more than 20 minutes of continuous use to cool down it, so that service life can be prolonged.
- The massage head is advised to be used in direct contact with skin surface. Healing effect could greatly diminish if used through clothing.
- When the massage head is pressed too hard on the skin, it will automatically stop vibrating to prevent muscle damage.

Biomedical-Combination of Treatment Technology[#]


Deep tissue stimulation

Using ergonomic percussion frequency to stimulate the musculoskeletal muscle and fascia. (Promote metabolism and relieve muscle stiffness)

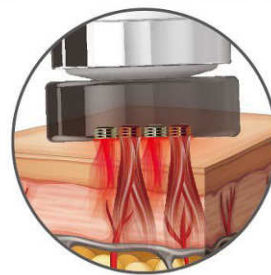Photobiomodulation

Research and development's low-level laser therapy system is composed of two medical wavelength signals and the energy level corresponds to medical standard. (Promote tissue repair, relieve pain and alleviate inflammation)

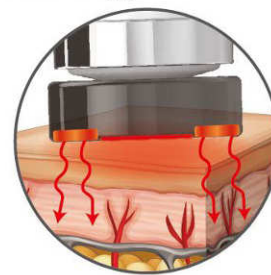FDA-reg heat therapy

Stimulation signal is controlled intelligently within FDA regulation on the temperature range (Promote rehabilitation and increase blood circulation)

Effectively relieve the following symptoms*

Sports injury
(Recommend to use after 24-36 hours)

Muscle strain

Arthritis/
Degenerative Diseases

- [#] Biomedical-Combination of Treatment Technology : A hybrid therapy bases on multiple therapeutic signals which the doses and formulations are exclusively designed by MTcure
- [#] Medical stimulation signals are designed according to medical standards with FDA's approval, please follow user manual's instructions strictly, especially the safety precautions' section 3 and 4.
- ^{*} Musculoskeletal conditions and related symptoms can be improved significantly from the clinical human test conducted by MTcure in collaboration with the Faculty of Medicine, The Chinese University of Hong Kong.